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NEBUR is an instrument suited for getting rid of syringe needles used to inject anaesthetic substances. It successfully eliminates any residues and avoid contamination thus ensuring to operate under the highest hygiene conditions. The product is available as a basic version or a **PLUS** version which offers the possibility to heat tubular vials before injection.

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
The appliance is mass-manufactured by C.I.E. DENTALFARM in compliance with technical and safety rules in force, as provided for by the Machinery Directive 89/392 EEC. Careful inspection and full routine testing is carried out singularly on each machine which is furtherly processed by an automatic testing installation (Model TS48 Elektrotech. Laboratorium, D-7015 KORNTAL 1) carrying out in sequence: Dielectric Strength Test, Test on the Protection Conductor, Insulation Resistance Test Short Circuit Test and it prints a report assuring compliance with the fixed limits.

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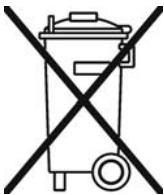According to International regulations, this unit has been classified as AEE (electric and electronic device, whose correct operation depends on electric currents and electromagnetic fields) and as a consequence, at the end of its lifetime, it cannot be treated as normal waste material but it must be disposed separately, complying with the Directive 2002/96/CE.

No special maintenance is required for correct use of the unit; only, periodically clean electrodes, visible in the opening, by means of a specific plastic tool. Elimination of the needle occurs by melting of the material it is made of, and fume emissions are thus present; at the same time, the liquid inside the needle will evaporate producing steam. Even if such vapours are in very small quantities, do not breath them and keep your face away from the hole while destroying the needle. To optimize battery life, read carefully the advice on this manual and do not continue using the battery when the charge level is low. Due to the presence of an overcharge control, the battery cannot be damaged by overcharging, but can suffer alteration when the battery level is low. Burnt residues will fall into the hole, but, even if they are sterile, we suggest to regularly empty the container by turning the cap leftwards, extract it and eliminate the residues as hospital and medical waste material. NOTE: To carry out this operation, do not capsize the unit and prevent residues from getting in touch with internal electronic components. No indications about battery life are available. One of the main symptoms of a used battery is its reduced capacity even after long charging cycles. NOTE: In case this inconvenience should repeat, check the unit with an authorized Technician.

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Insert the Anesthetic vial into the designated opening, press the Red function key. The red Led will light up initiating the warm-up for 8-9 minutes heating stage. After 2 minutes the Led will switch off, indicating that the temperature has reached 38-39°C. At this stage, the anesthetic is warm and ready for use; since that moment, the warm up chamber will keep this temperature for 8-9 minutes (red LED goes on and off regularly). To prolong the heating stage, press the Red key again to set additional 8-9 minute heating time. For battery saving purposes, do not use this function too frequently.

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Read carefully the present manual before carrying out any starting, operating or maintenance operation on the machine

Safety of the User depends on the skill ability, on the good sense and on the caution when using the machine; for this reason, it is of utmost importance to know in detail the allocation and the function of all the controls.

The final User of the machine is a skilled and well trained professional and not merely an occasional worker

Check regularly the parts which tend to wear out due to the specific working conditions

Do not tamper the electrical wiring system of the machine

The built-in safety mechanisms shall in no way be removed or modified

Do not allow unauthorised persons to try to repair the machine

Do not utilise petrol, solvents or any other inflammable fluids to clean some parts of the machine – use exclusively non-toxic and non-inflammable products available on the market for this specific purpose.

Dentalfarm is no way responsible in case above safety precautions are not strictly observed.

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Press button «test batt.» above the unit, make sure the nearby colourless Led lights up in YELLOW/GREEN (the same happens when connected to the Battery Charger); this means the battery is sufficiently charged for use. If the same Led lights up Red, do NOT use the device, as this states for too low battery for use.

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Make sure the voltage indicated in this booklet as well as on the Battery Charger power plate corresponds to your mains voltage.

Connect the power jack of the Battery Charger to a power socket; make sure the power outlet is properly earthed in compliance with applicable regulations.

Fit the low voltage jack into the power socket under the device (to do this, do not capsize the device since the combustion residues in the container would fall out; as a consequence, we suggest to empty the container every time the battery is charged, before fitting the jack).

The Green «charge batt» Led will switch on and will keep on indicating the device is connected to the power mains.

This condition is not a hazard warning for the operator and it does not damage the device. The unit will continue operating perfectly under this condition.

For correct battery use and in order to ensure long and efficient battery life, we recommend not to use the device permanently under this condition.

After about 12 hours, the colourless Led in off position will switch on (YELLOW/GREEN light) to indicate that the battery is fully charged.

Remove the Battery Charger from the device and the mains. The device features a maximum battery charge control so that even if this operation is not done in the immediate, neither the battery nor the device can be damaged.

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The operating principle is the same of the electric arc; the needle inserted in the hole closes the electrical contact between the two battery poles, whose electric shock is higher than that of the needle electrical resistance.

The needle immediately heats up, over 1,400°C and melts and it is thus eliminated. The remaining part, besides being thermally sterilized, is no longer sharp and is melt, preventing any leakage of liquid.

When the battery level is high enough, the device is constantly ready for use and, before next charge, can successfully eliminate around 200/800 needles depending on their size.

NEBUR NEEDLE BURNER is able to destroy:

- Needles of any kind: intramuscular, insuline syringes, anaesthesia and endocanal injections, up to 15 Gauge.

- Endodontic instruments: files, reamers.

- Small cutting instruments: scalpel blades and all small pointed or sharp disposable instruments.

WORKING RECOMMENDATIONS

Immediately after use: whether the needle has to be removed or not, no matter which instrument is used (Carpule, Peripress or Cito Ject syringes or other supports), before removing and/or discarding the object, insert it into the NEBUR opening in order to exclude any possibility of infection. This operation takes up to 2-3 seconds in average.